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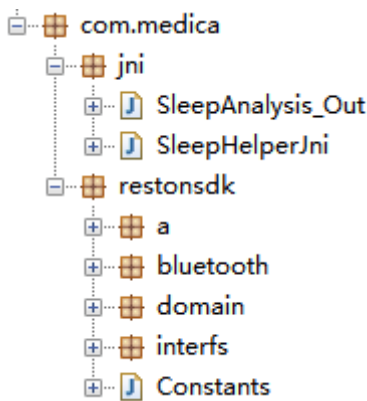
public short MdRemSleepPerc
中睡眠百分

- **MdLightSleepPerc**

public short MdLightSleepPerc
浅睡眠百分

- **MdWakeSleepPerc**

public short MdWakeSleepPerc
清醒睡眠百分比

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public short MdHeartRateMin
最低心率

- **MdBreathRateMin**

public short MdBreathRateMin
最低呼吸率

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public short md_body_move_decrease_scale
躁动不安扣分

- **md_leave_bed_decrease_scale**

public short md_leave_bed_decrease_scale

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- **md_wake_cnt_decrease_scale**

public short md_wake_cnt_decrease_scale

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- **md_start_time_decrease_scale**

public short md_start_time_decrease_scale

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- **md_fall_asleep_time_decrease_scale**

public short md_fall_asleep_time_decrease_scale

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- **md_perc_deep_decrease_scale**

public short md_perc_deep_decrease_scale

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- **md_is_sleep_time_long**

public short md_is_sleep_time_long

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- **md_sleep_time_decrease_scale**

public short md_sleep_time_decrease_scale

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- **md_breath_stop_decrease_scale**

public short md_breath_stop_decrease_scale

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- **md_heart_stop_decrease_scale**

public short md_heart_stop_decrease_scale

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- **md_heart_low_decrease_scale**

public short md_heart_low_decrease_scale

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- **md_heart_high_decrease_scale**

public short md_heart_high_decrease_scale

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- **md_breath_low_decrease_scale**

public short md_breath_low_decrease_scale

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- **md_breath_high_decrease_scale**

public short md_breath_high_decrease_scale

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- **STATE_DISCONNECTED**

```
public static final int STATE_DISCONNECTED  
ble 设备连接状态: 连接断开
```

- **STATE_CONNECTING**

```
public static final int STATE_CONNECTING  
ble 设备连接状态: 正在连接
```

- **STATE_CONNECTED**

```
public static final int STATE_CONNECTED  
ble 设备连接状态: 已连接
```

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- **onStateChanged**

void onStateChanged(int state)
参数：
state - 取值参考
RestOnHelper.STATE_DISCONNECTED,
RestOnHelper.STATE_CONNECTING,
RestOnHelper.STATE_CONNECTED

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- **CONNECT_DEVICE**

public static final Method CONNECT_DEVICE

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- **LOGIN**

public static final Method LOGIN

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- **LOGOUT**

public static final Method LOGOUT

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- **GET_DEVICE_VERSION**

public static final Method GET_DEVICE_VERSION
获取设备版本

- **GET_DEVICE_STATUS**

public static final Method GET_DEVICE_STATUS
获取设备状态

- **START_COLLECT**

public static final Method START_COLLECT
开始采集数据

- **STOP_COLLECT**

public static final Method STOP_COLLECT
停止采集数据

- **SEE_REALTIME_DATA**

public static final Method SEE_REALTIME_DATA
查看实时数据

- **CLOSE_REALTIME_DATA**

public static final Method CLOSE_REALTIME_DATA
关闭实时数据

- **SEE_RAW_DATA**

public static final Method SEE_RAW_DATA
查看原始数据

- **CLOSE_RAW_DATA**

public static final Method CLOSE_RAW_DATA
关闭原始数据

- **QUERY_HISTORY_SUMMARY**

public static final Method QUERY_HISTORY_SUMMARY
查询历史数据概要信息

- **QUERY_HISTORY_DETAIL**

public static final Method QUERY_HISTORY_DETAIL
查询历史数据详细信息

- **GET_DEVICE_POWER**

public static final Method GET_DEVICE_POWER
获取 RestOn 电量百分比

- **CLOSE_LOWPOWER_WARNING**

public static final Method CLOSE_LOWPOWER_WARNING
关闭低电量警告

- **SET_AUTO_START**

public static final Method SET_AUTO_START
设置自动开始监测

- **SET_ALARM_TIME**

public static final Method SET_ALARM_TIME
设置闹钟时间

- **GET_DEVICE_ID**

public static final Method GET_DEVICE_ID
获取设备 id

- **UPGRADE_FIRMWARE**

public static final Method UPGRADE_FIRMWARE
升级固件

```
package com.medica.restonsdk;  
/**  
 * RestOn SDK相关常量  
 */  
public class Constants
```

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- **SDK_VERSION**

public static final String SDK_VERSION
SDK 版本

- **CONFIG_FILENAME**

public static final String CONFIG_FILENAME
配置文件名称

- **DEVICE_TYPE_RESTON_Z1**

public static final short DEVICE_TYPE_RESTON_Z1
设备型号 RestOn Z1

- **DEVICE_TYPE_RESTON_Z2**

```
public static final short DEVICE_TYPE_RESTON_Z2
设备型号 RestOn Z2
```

```
package com.medica.restonsdk;
/**
 * 监测状态值
 */
public class Constants.SleepStatusType
```

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```
public static final byte SLEEP_BODYMOVE
体动
```

• SLEEP_LEAVE

```
public static final byte SLEEP_LEAVE
离床
```

• SLEEP_TURN_OVER

```
public static final byte SLEEP_TURN_OVER
翻身
```

• SLEEP_BODYMOVE_TEMP

```
public static final byte SLEEP_BODYMOVE_TEMP
体动振幅
```

- SLEEP_INVALID

public static final byte SLEEP_INVALID

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com.medica.restonsdk.Constants		
限定符和类型	常量字段	值
public static final <u>String</u>	<u>CONFIG_FILENAME</u>	"restonSdkConfig"
public static final short	<u>DEVICE_TYPE_RESTON_Z1</u>	1
public static final short	<u>DEVICE_TYPE_RESTON_Z2</u>	9
public static final <u>String</u>	<u>SDK_VERSION</u>	"1.0.0"

com.medica.restonsdk.Constants.SleepStatusType		
限定符和类型	常量字段	值
public static final byte	<u>SLEEP_B_STOP</u>	2
public static final byte	<u>SLEEP_BODYMOVE</u>	4
public static final byte	<u>SLEEP_BODYMOVE_TEMP</u>	7
public static final byte	<u>SLEEP_H_STOP</u>	3
public static final byte	<u>SLEEP_INIT</u>	1
public static final byte	<u>SLEEP_INVALID</u>	-1
public static final byte	<u>SLEEP_LEAVE</u>	5
public static final byte	<u>SLEEP_OK</u>	0
public static final byte	<u>SLEEP_TURN_OVER</u>	6

com.medica.restonsdk.bluetooth.RestOnHelper		
限定符和类型	常量字段	值
public static final <u>String</u>	<u>ACTION_LOW_POWER</u>	"com.medica.restonsdk.Device_LowPower"
public static final int	<u>STATE_CONNECTED</u>	2
public static final int	<u>STATE_CONNECTING</u>	1
public static final int	<u>STATE_DISCONNECTED</u>	0